


Penguin Switch Adapted Toy

USER GUIDE

Introduction

The Penguin Switch Adapted Toy sings and dances and can be controlled with accessible switches that have a 3.5 mm jack. The switch adapted toy can be used to get the penguin to do the alphabet, opposites, counting, and bob side to side with lights and music.

Features

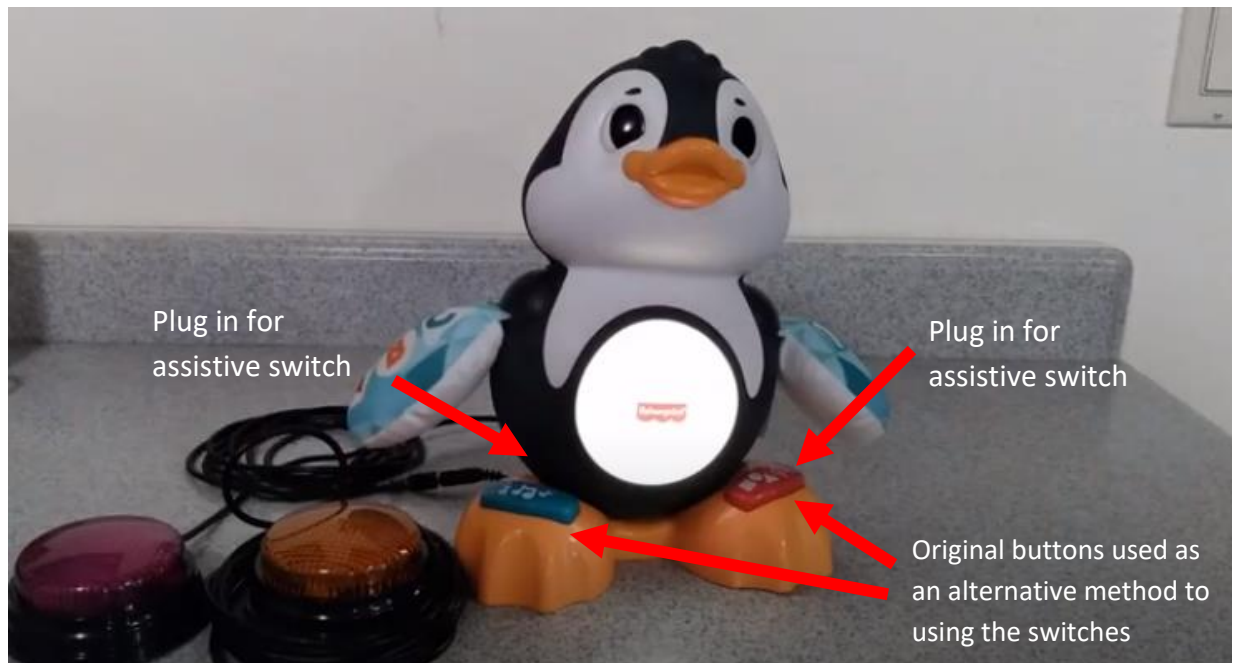

Usage

To activate the penguin, the power switch must be turned on, which is located on bottom of the penguin. Once turned on, the assistive switches can be used as an alternative to the green and red buttons on the penguin's feet for activation. The switch adapted toy would be beneficial for someone who may have difficulty pressing the original buttons on the toy. Some people may find the original buttons too hard to press, too small, or out of reach for appropriate use.

Compatibility

The Penguin Switch Adapted Toy is compatible with any switches that have 3.5 mm cable ends.

Specifications

Height (mm)	Width (mm)	Depth (mm)	Weight (g)
270	247	138	970

Cleaning

Wipe the toy with a clean cloth. Be careful around the exposed electronic jacks if the cloth is damp.